

### Preset Scenarios:

- F1 – SNNN
- F2 – SSNN
- F3 – NSNN
- F4 – NNNN
- S=Shock advised
- N=No shock advised
- C1 – Custom #1: SSSSSSSS
- C2 – Custom #2: SSSSSSSS


### Other options:

- OP – Sets operating settings for trainer
- T1 – CPR time 1
- T2 – CPR time 2
- L – Language

### During Training Session:

- Not used for CR-T
- Simulates pads are placed correctly and will move on to analyzing the patient's heart rhythm
- Changes "SHOCK ADVISED" to "NO SHOCK ADVISED" during analysis
- Displays CHARGE-PAK symbol to simulate need to replace CHARGE-PAK™ battery charger
- Simulates low internal battery condition
- Simulates a need for service
- Simulates motion and will prompt the user to stop moving the patient during the analysis
- Pause
- Adjusts volume

Please refer to the CR-T Trainer Operating Instructions for complete directions for use.

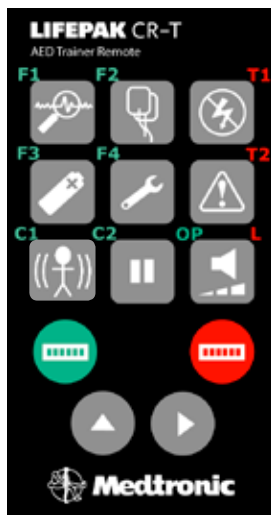

## LIFEPAK® CR-T AED TRAINING SYSTEM

### Setting Up Your CR-T



**Physio-Control Headquarters**  
 11811 Willows Road NE  
 Redmond, WA 98052  
[www.physio-control.com](http://www.physio-control.com)

*Customer Support*  
 P. O. Box 97006  
 Redmond, WA 98073  
 Toll Free 800 442 1142  
 Fax 800 426 8049

**Physio-Control Canada**  
 Medtronic of Canada Ltd  
 99 Hereford Street  
 Brampton, ON  
 L6Y 0R3  
 Tel 888 879 0977  
 Fax 866 430 6115

**Physio-Control, Inc.**, 11811 Willows Road NE,  
 Redmond, WA 98052 USA






## Setting Up Your CR-T

Slide the power switch on the right side of the device forward to turn on device. OK symbol should now be in the readiness display.





Use the remote control to set your device into the following settings:

### Sequence - Shock then No Shock Advised (Fixed setting #1)

- Press  on the remote
- Press **F1** (  )
- On the LED display you should see the first bars as green and the second, third and fourth bars as red
- Press  on the remote. It is now set for Shock, No Shock Advised.





To configure other shock scenarios, go to Page 19 in the *CR-T* Operating Instructions.

### Set to Fully-Automatic (no use of shock button) mode

- Press  on the remote
- Press **OP** (  )
- Use  to make the first bar in the LED display red
- Press 
- Be sure to put plastic cover over shock button on device

To configure for the semi-automatic setting, go to page 20 in the *CR-T* Operating Instructions.

### For a faster demo, set CPR time to 15 seconds

- Press  on the remote
- Press **T1** (  )
- Move the red bar in the LED display to the first position using the  button
- Press  on the remote

To set other CPR time options, go to page 21 in the *CR-T* Operating Instructions.

### Tip:

When using the remote, be sure it is within range of the LED display/IRDA port so it will read commands.